# Operating Instructions **Memosens CLS15E**

Conductivity sensor with Memosens protocol For conductive measurement of conductivity in liquids



BA02018C/07/EN/01.20

71515065 2021-01-31





# Table of contents

1	About this document 3
1.1	Warnings 3
1.2	Symbols 3
1.3	Documentation 3
2	Basic safety instructions 4
2.1	Requirements for personnel 4
2.2	Designated use 4
2.3	Workplace safety 4
2.4	Operational safety 5
2.5	Product safety 5
3	Incoming acceptance and
	product identification 5
3.1	Incoming acceptance 5
3.2	Product identification 6
3.3	Scope of delivery 6
3.4	Certificates and approvals 7
4	Installation7
4.1	Installation conditions 7
4.2	Mounting the sensor
4.3	Post-installation check
5	Electrical connection
5.1	Connecting the sensor 10
5.2	Ensuring the degree of protection 10
5.3	Post-connection check 10
6	Commissioning 11
7	Maintenance 11
8	Repair 12
8.1	General information
8.2	Spare parts
8.3	Return
8.4	Disposal 13
9	Accessories 14
9.1	Threaded and adapter couplings 14
9.2	Flow assembly
9.3	Measuring cable 14
9.4	Calibration solutions 15
9.5	Calibration set 15

Memosens CL	S15E
-------------	------

10	Technical data 15	
10.1	Input 15	
10.2	Performance characteristics 16	
10.3	Environment 16	
10.4	Process 17	
10.5	Mechanical construction 17	
Index 19		

### 1 About this document

### 1.1 Warnings

Structure of information	Meaning
▲ DANGER Causes (/consequences) If necessary, Consequences of non- compliance (if applicable) ► Corrective action	This symbol alerts you to a dangerous situation. Failure to avoid the dangerous situation <b>will</b> result in a fatal or serious injury.
WARNING Causes (/consequences) If necessary, Consequences of non- compliance (if applicable) ► Corrective action	This symbol alerts you to a dangerous situation. Failure to avoid the dangerous situation <b>can</b> result in a fatal or serious injury.
▲ CAUTION Causes (/consequences) If necessary, Consequences of non- compliance (if applicable) ► Corrective action	This symbol alerts you to a dangerous situation. Failure to avoid this situation can result in minor or more serious injuries.
NOTICE Cause/situation If necessary, Consequences of non- compliance (if applicable) Action/note	This symbol alerts you to situations which may result in damage to property.

### 1.2 Symbols

Symbol	Meaning
i	Additional information, tips
	Permitted or recommended
	Not permitted or not recommended
I	Reference to device documentation
	Reference to page
	Reference to graphic
۲	Result of a step

### 1.3 Documentation

The following manual which complements these Operating Instructions can be found on the product page on the Internet:

Technical Information Memosens CLS15E, TI01526C

In addition to these Operating Instructions, an XA with "Safety instructions for electrical apparatus in the hazardous area" is also included with sensors for use in the hazardous area.

▶ Please follow instructions on use in the hazardous area carefully.

## 2 Basic safety instructions

### 2.1 Requirements for personnel

- Installation, commissioning, operation and maintenance of the measuring system may be carried out only by specially trained technical personnel.
- The technical personnel must be authorized by the plant operator to carry out the specified activities.
- The electrical connection may be performed only by an electrical technician.
- The technical personnel must have read and understood these Operating Instructions and must follow the instructions contained therein.
- Faults at the measuring point may only be rectified by authorized and specially trained personnel.

Repairs not described in the Operating Instructions provided must be carried out only directly at the manufacturer's site or by the service organization.

### 2.2 Designated use

The conductivity sensor is designed for the conductive measurement of the conductivity of liquids.

It is used in the following areas:

Measurements in pure and ultrapure water

Use of the device for any purpose other than that described, poses a threat to the safety of people and of the entire measuring system and is therefore not permitted.

The manufacturer is not liable for damage caused by improper or non-designated use.

### 2.3 Workplace safety

As the user, you are responsible for complying with the following safety conditions:

- Installation guidelines
- Local standards and regulations
- Regulations for explosion protection

#### Electromagnetic compatibility

- The product has been tested for electromagnetic compatibility in accordance with the applicable international standards for industrial applications.
- The electromagnetic compatibility indicated applies only to a product that has been connected in accordance with these Operating Instructions.

### 2.4 Operational safety

#### Before commissioning the entire measuring point:

- 1. Verify that all connections are correct.
- 2. Ensure that electrical cables and hose connections are undamaged.
- 3. Do not operate damaged products, and protect them against unintentional operation.
- 4. Label damaged products as defective.

#### During operation:

 If faults cannot be rectified: products must be taken out of service and protected against unintentional operation.

#### 2.5 Product safety

The product is designed to meet state-of-the-art safety requirements, has been tested, and left the factory in a condition in which it is safe to operate. The relevant regulations and international standards have been observed.

### 3 Incoming acceptance and product identification

#### 3.1 Incoming acceptance

- 1. Verify that the packaging is undamaged.
  - Notify the supplier of any damage to the packaging.
    Keep the damaged packaging until the issue has been resolved.
- 2. Verify that the contents are undamaged.
  - Notify the supplier of any damage to the delivery contents.
    Keep the damaged goods until the issue has been resolved.
- 3. Check that the delivery is complete and nothing is missing.
  - └ Compare the shipping documents with your order.
- 4. Pack the product for storage and transportation in such a way that it is protected against impact and moisture.
  - The original packaging offers the best protection.
    Make sure to comply with the permitted ambient conditions.

If you have any questions, please contact your supplier or your local Sales Center.

### 3.2 Product identification

#### 3.2.1 Nameplate

The nameplate provides you with the following information on your device:

- Manufacturer identification
- Extended order code
- Serial number
- Safety information and warnings
- Compare the information on the nameplate with the order.

#### 3.2.2 Identifying the product

#### Product page

www.endress.com/cls15e

#### Interpreting the order code

The order code and serial number of your product can be found in the following locations:

- On the nameplate
- In the delivery papers

#### Obtaining information on the product

- 1. Go to www.endress.com.
- 2. Call up the site search (magnifying glass).
- 3. Enter a valid serial number.
- 4. Search.
  - └ The product structure is displayed in a popup window.
- 5. Click on the product image in the popup window.
  - ← A new window (**Device Viewer**) opens. All of the information relating to your device is displayed in this window as well as the product documentation.

#### Manufacturer address

Endress+Hauser Conducta GmbH+Co. KG Dieselstraße 24 D-70839 Gerlingen

### 3.3 Scope of delivery

The scope of delivery includes:

- Sensor in the version ordered
- Operating Instructions

### 3.4 Certificates and approvals

Certificates and approvals are optional, i.e. they depend on the product version.

#### 3.4.1 C€ mark

#### EU Declaration of Conformity

The product meets the requirements of the harmonized European standards. As such, it complies with the legal specifications of the EU directives. The manufacturer confirms successful testing of the product by affixing to it the CC mark.

#### 3.4.2 Test reports

#### Manufacturer inspection certificate

Stating the individual cell constant

#### 3.4.3 Additional certification

#### Inspection certificate in accordance with EN 10204 3.1

A test certificate 3.1 in accordance with EN 10204 is supplied depending on the version ( $\rightarrow$  Product Configurator on the product page).

#### 3.4.4 Other standards and guidelines

#### EAC


The product has been certified according to guidelines TP TC 004/2011 and TP TC 020/2011 which apply in the European Economic Area (EEA). The EAC conformity mark is affixed to the product.

## 4 Installation


### 4.1 Installation conditions

The sensors are installed directly via the process connection, e.g. Clamp.

As an option, the sensor can also be installed using a commercially available T-piece or cross fitting or using a flow assembly.

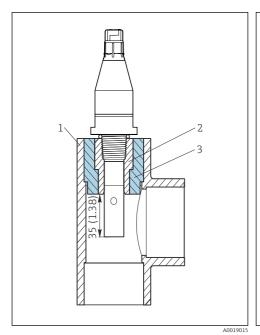


- With NPT ½" thread in T-piece or cross fitting
- 1 T-piece or cross fitting (DN 32, 40 or 50)
- 2 Glue-in VC threaded coupling (NPT ½" for DN 20)
- 3 Glue-in adapter coupling (for DN 32, 40, 50)

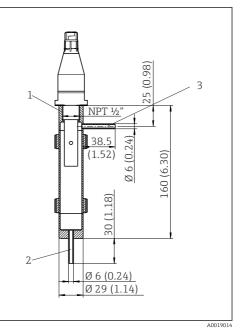

- With NPT ½" thread in flow assembly 71042405, dimensions in mm (inch)
  - Sensor holder NPT ½"
- 2 Inlet
- 3 Outlet
- ► If using the sensor in the ultrapure water range, you must work under air-evacuated conditions.
  - $\vdash$  Otherwise, the CO<sub>2</sub> in the air can dissolve in the water and its (weak) dissociation can increase the conductivity by up to 3 µS/cm.

1

#### 4.2 Mounting the sensor

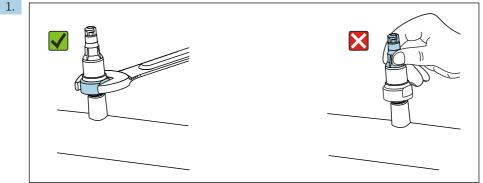

A004290

#### NOTICE

#### Incorrect mounting or disassembly

The head could become loose and fall off, resulting in total sensor failure!

- Only mount the sensor via the process connection.
- ► To do so, use a suitable tool, such as an open-ended wrench.

Install the sensor via the process connection or an assembly.

2. Ensure that the electrodes are fully immersed in the medium during measurement. Immersion depth: at least 35 mm (1.38").

### 4.3 Post-installation check

- 1. Are the sensor and cable undamaged?
- 2. Is the sensor installed in the process connection and is not suspended from the cable?

# 5 Electrical connection

#### **WARNING**


#### Device is live!

Incorrect connection may result in injury or death!

- ► The electrical connection may be performed only by an electrical technician.
- The electrical technician must have read and understood these Operating Instructions and must follow the instructions contained therein.
- ▶ **Prior** to commencing connection work, ensure that no voltage is present on any cable.

### 5.1 Connecting the sensor

The electrical connection of the sensor to the transmitter is established using the measuring cable CYK10.

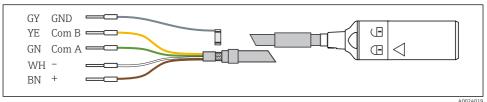

☑ 3 Measuring cable CYK10

#### NOTICE

#### Mechanical twist protection

If too much force is applied to the Memosens head, this can shear the connections and therefore destroy the sensor!

- ► There is no need to exert excessive force when connecting the sensor to the cable coupling. Proceed with care!
- ► If the Memosens coupling clearly will not close, check the coupling for dirt or mechanical damage, and make sure that you are turning it in the right direction. Pay attention to the lock symbol on the coupling!
- ▶ If necessary, use a different Memosens cable.

### 5.2 Ensuring the degree of protection

Only the mechanical and electrical connections which are described in these instructions and which are necessary for the required, designated use, may be carried out on the device delivered.

• Exercise care when carrying out the work.

Otherwise, the individual types of protection (Ingress Protection (IP), electrical safety, EMC interference immunity) agreed for this product can no longer be guaranteed due, for example to covers being left off or cable (ends) that are loose or insufficiently secured.

### 5.3 Post-connection check

#### **WARNING**

#### **Connection errors**

The safety of people and of the measuring point is at risk! The manufacturer does not accept any responsibility for errors that result from failure to comply with the instructions in this manual.

 Put the measuring point into operation only if you can answer yes to all the following questions.

Product status and specifications

• Are the sensor and cable free from damage on the outside?

Electrical connection

- ▶ Is the installed cable strain-relieved and not twisted?
- Is a sufficient length of the cable cores stripped, and are the cores correctly positioned in the terminal on the transmitter?
- ► Are all plug-in terminals on the transmitter securely engaged?
- ▶ Are all cable entries mounted on the transmitter, tightened and leak-tight?

#### 6 Commissioning

Prior to initial commissioning, ensure that:

- The sensor is correctly installed
- The electrical connection is correct



1. Check the temperature compensation and damping settings on the transmitter.



Operating Instructions of the transmitter used, e.g. BA01245C if the Liquiline CM44x or CM44xR is used.

#### **WARNING**

#### Escaping process medium

Risk of injury from high pressure, high temperatures or chemical hazards!

- Before applying pressure to an assembly with cleaning system, ensure that the system has been connected correctly.
- ▶ If you cannot reliably establish the correct connection, do not install the assembly in the process.

If using an assembly with automatic cleaning function:

- 2. Check that the cleaning medium (water or air, for example) is connected correctly.
- 3. Following commissioning:

Maintain the sensor at regular intervals.

└ This is the only way to ensure a reliable measurement.

#### 7 Maintenance

#### **A**CAUTION

#### **Corrosive chemicals**

Risk of chemical burns to the eyes and skin and risk of damage to clothing and equipment!

- ▶ It is absolutely essential to protect the eyes and hands properly when working with acids, alkalis and organic solvents!
- Wear protective goggles and safety gloves.
- ► Clean away splashes on clothes and other objects to prevent any damage.
- ► Comply with instructions in the safety data sheets for the chemicals used.

### **WARNING**

#### Thiocarbamide

Harmful if swallowed! Limited evidence of carcinogenicity! Possible risk of harm to the unborn child! Dangerous for the environment with long-term effects!

- ► Wear protective goggles, protective gloves and appropriate protective clothing.
- Avoid all contact with the eyes, mouth and skin.
- Avoid discharge into the environment.

Clean away fouling on the sensor as follows depending on the type of fouling:

1. Oily and greasy films:

Clean with fat solvent, e.g. alcohol, or hot water and agents containing surfactants (alkaline) (e.g. dishwashing detergent).

- 2. Lime and metal hydroxide buildup and low solubility (lyophobic) organic buildup: Dissolve buildup with diluted hydrochloric acid (3 %) and then rinse thoroughly with plenty of clear water.
- Sulfidic buildup (from flue gas desulfurization or wastewater treatment plants):
  Use a mixture of hydrochloric acid (3 %) and thiocarbamide (commercially available) and then rinse thoroughly with plenty of clear water.
- Buildup containing proteins (e.g. food industry):
  Use a mixture of hydrochloric acid (0.5 %) and pepsin (commercially available) and then rinse thoroughly with plenty of clear water.
- 5. Readily soluble biological buildup: Rinse with pressurized water.

After cleaning, rinse the sensor thoroughly with water.

# 8 Repair

### 8.1 General information

The repair and conversion concept provides for the following:

- The product has a modular design
- Spare parts are grouped into kits which include the associated kit instructions
- Only use original spare parts from the manufacturer
- Repairs are carried out by the manufacturer's Service Department or by trained users
- Certified devices can only be converted to other certified device versions by the manufacturer's Service Department or at the factory
- Observe applicable standards, national regulations, Ex documentation (XA) and certificates
- 1. Carry out the repair according to the kit instructions.
- 2. Document the repair and conversion and enter, or have entered, in the Life Cycle Management tool (W@M).

### 8.2 Spare parts

Device spare parts that are currently available for delivery can be found on the website:

www.endress.com/device-viewer

• Quote the serial number of the device when ordering spare parts.

#### 8.3 Return

The product must be returned if repairs or a factory calibration are required, or if the wrong product was ordered or delivered. As an ISO-certified company and also due to legal regulations, Endress+Hauser is obliged to follow certain procedures when handling any returned products that have been in contact with medium.

To ensure the swift, safe and professional return of the device:

 Refer to the website www.endress.com/support/return-material for information on the procedure and conditions for returning devices.

#### 8.4 Disposal



If required by the Directive 2012/19/EU on waste electrical and electronic equipment (WEEE), the product is marked with the depicted symbol in order to minimize the disposal of WEEE as unsorted municipal waste. Do not dispose of products bearing this marking as unsorted municipal waste. Instead, return them to Endress+Hauser for disposal under the applicable conditions.

### 9 Accessories

The following are the most important accessories available at the time this documentation was issued.

► For accessories not listed here, please contact your Service or Sales Center.

### 9.1 Threaded and adapter couplings

For sensors with NPT <sup>1</sup>/<sub>2</sub>" process connection

#### PVC threaded coupling

- For gluing into commercially available PVC cross fittings or T-pieces with DN 20
- With G<sup>1</sup>/<sub>2</sub> internal thread, self-sealing with NPT <sup>1</sup>/<sub>2</sub>" sensor thread
- Order No. 50066536

#### PVDF threaded coupling

- With G<sup>1</sup>/<sub>2</sub> internal thread and G1 external thread
- Explosion-proof up to 12 bar at 20 °C (174 psi at 68 °F), max. 120 °C at 1 bar (248 °F at 14.5 psi), including O-ring
- Self-sealing internal thread with NPT <sup>1</sup>/<sub>2</sub>" sensor thread
- Order No. 50004381

#### PVC adapter couplings AM

- For adapting the PVC threaded coupling to larger nominal diameters
- Diameter, order numbers:
  - AM 32: for cross-fittings or T-pieces DN 32, Order No. 50004738
  - AM 40: for cross-fittings or T-pieces DN 40, Order No. 50004739
  - AM 50: for cross-fittings or T-pieces DN 50, Order No. 50004740

### 9.2 Flow assembly

For installing conductivity sensors with NPT 1/2" thread

- Material: Stainless steel 1.4404 (AISI 316 L)
- Connections: 90°, Ø 6 mm (0.24")
- Volume: 0.69 l (0.18 US gal)
- Max. temperature: 100 °C (212 °F)
- Max. pressure: 16 bar (232 psi)
- Order No.: 71042405

### 9.3 Measuring cable

#### Memosens data cable CYK10

- For digital sensors with Memosens technology
- Product Configurator on the product page: www.endress.com/cyk10

Technical Information TI00118C

#### Memosens data cable CYK11

- Extension cable for digital sensors with Memosens protocol
- Product Configurator on the product page: www.endress.com/cyk11



Technical Information TI00118C

#### 9.4 **Calibration solutions**

#### Conductivity calibration solutions CLY11

Precision solutions referenced to SRM (Standard Reference Material) by NIST for qualified calibration of conductivity measuring systems in accordance with ISO 9000

- CLY11-A, 74 μS/cm (reference temperature 25 °C (77 °F)), 500 ml (16.9 fl.oz) Order No. 50081902
- CLY11-B, 149.6 μS/cm (reference temperature 25 °C (77 °F)), 500 ml (16.9 fl.oz) Order No. 50081903



Technical Information TI00162C

#### 9.5 **Calibration** set

#### Conducal CLY421

- Conductivity calibration set (case) for ultrapure water applications
- Complete, factory-calibrated measuring system with certificate, traceable to SRM by NIST and PTB, for comparison measurement in ultrapure water up to max. 20 µS/cm
- Product Configurator on the product page: www.endress.com/cly421

Technical Information TI00496C/07/EN

#### Technical data 10

#### 10.1 Input

#### 10.1.1 Measured variables

- Conductivity
- Temperature

#### Measuring ranges 10.1.2

#### Conductivity<sup>1)</sup>

CLS15E-\*\*\*\*A CLS15E-\*\*\*\*B

40 nS/cm to 20 µS/cm 100 nS/cm to 200  $\mu$ S/cm

1) In relation to water at 25 °C (77 °F)

#### Temperature

-20 to 140 °C (-4 to 284 °F)

#### 10.1.3 Cell constant

CLS15E-*****A	$k = 0.01 \text{ cm}^{-1}$
CLS15E-*****B	$k = 0.1 \text{ cm}^{-1}$

#### 10.1.4 Temperature compensation

Pt1000 (Class A according to IEC 60751)

#### 10.2 Performance characteristics

#### 10.2.1 Uncertainty of measurement

Each individual sensor is factory-measured in a solution with approx. 5  $\mu$ S/cm for cell constant 0.01 cm<sup>-1</sup> or approx. 50  $\mu$ S/cm for cell constant 0.1 cm<sup>-1</sup> using a reference measuring system traceable to NIST or PTB. The exact cell constant is entered into the manufacturer inspection certificate supplied. The uncertainty of measurement in determining the cell constant is 1.0 %.

#### 10.2.2 Response time

Conductivity	t <sub>95</sub> ≤ 2 s
Temperature <sup>1)</sup>	
CLS15E-*****A:	$t_{90} \le 16 \text{ s}^{2)}$
CLS15E-****B:	$t_{90} \le 8 \ s^{2}$

1) DIN VDI/VDE 3522-2 (0.3 m/s laminar)

2) With temperature prediction activated as standard

#### 10.2.3 Measured error

Conductivity	$\leq$ 2 % of reading, in specified measuring range
-	$\leq$ 0.5 K, in measuring range -5 to 100 °C (23 to 212 °F) $\leq$ 1.0 K, in measuring range 100 to 140 °C (212 to 284 °F)

#### 10.2.4 Repeatability

Conductivity	$\leq 0.2$ % of reading, in specified measuring range
Temperature	≤ 0.05 K

#### 10.3 Environment

#### 10.3.1 Ambient temperature

-20 to 60 °C (-4 to 140 °F)

#### 10.3.2 Storage temperature

-25 to +80 °C (-10 to +180 °F)

#### 10.3.3 Degree of protection

IP 68 / NEMA type 6P (1.9 m water column, 20 °C, 24 h)

#### 10.4 Process

#### 10.4.1 Process temperature


Normal operation	-20 to 120 °C (-4 to 248 °F)
Sterilization (max. 1 h) $^{1)}$	Max. 140 °C (284 °F)

1) Threaded versions: max. 30 minutes

#### 10.4.2 Process pressure

13 bar (188 psi) absolute, at 20  $^{\circ}$ C (68  $^{\circ}$ F) 2 bar (29 psi) absolute, at 120  $^{\circ}$ C (248  $^{\circ}$ F)

#### 10.4.3 Temperature/pressure ratings

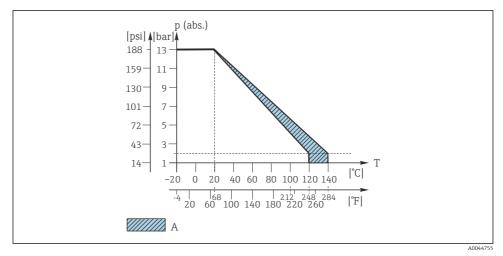

Mechanical pressure-temperature resistance

A Can be sterilized for a short time (1 hour)

#### 10.5 Mechanical construction

#### 10.5.1 Weight

Approx. 0.3 kg (0.66 lbs) depending on version

#### 10.5.2 Materials (in contact with medium)

Electrodes	Polished, stainless steel 1.4435 (AISI 316L)
Sensor shaft	Polyethersulfone (PES-GF20)
O-ring, in contact with medium	EPDM
(only Clamp version)	

#### 10.5.3 Materials (not in contact with medium)

#### Information according to REACH Regulation (EC) 1907/2006 Art. 33/1)

An internal connector contains the SVHC substance lead (CAS number 7439-92-1) with over 0.1 % (w/w).

The product does not present a hazard if it is used as designated.

#### 10.5.4 Process connection

Thread NPT  $\frac{1}{2}$ " and  $\frac{3}{4}$ " Clamp 1 $\frac{1}{2}$ " as per ISO 2852

#### 10.5.5 Surface roughness

 $R_a \le 0.8 \ \mu m$ 

# Index

# Α

Ambient temperature	16
С	
Cell constant	16
Check	
Connection	10
Installation	9
Connection	
Check	10
Ensuring the degree of protection	10

### D

Declaration of Conformity 7
Degree of protection
Ensuring
Technical data
Designated use
Disposal
Documentation

### E

Electrical connection	9
EN 10204 3.1	7
Environment	16

# I

Identifying the product	6
Incoming acceptance	5
Inspection certificate	7
Installation	
Check	9
Sensor	9

### М

Manufacturer inspection certificate	. 7
Materials	18
Measured error	16
Measured variables	15
Measuring ranges	15

### N

Nameplate				•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•	•		(	5
-----------	--	--	--	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	--	---	---

### 0

Operational safety												•					5	
--------------------	--	--	--	--	--	--	--	--	--	--	--	---	--	--	--	--	---	--

### Ρ

Performance characteristics	16
Pressure/temperature ratings	17
Process	17
Process connection	18
Process pressure	17
Process temperature	17
Product safety	5

### R

Repair	12
Repeatability	16
Response time	16
Return	13

# Safe

Sarety
Operation
Product
Workplace safety
Safety instructions
Scope of delivery
Sensor
Cleaning
Connecting
Mounting
Spare parts
Storage temperature
Surface roughness
Symbols

### Т

Technical data	
Environment	16
Input	15
Mechanical construction	17
Performance characteristics	16
Process	17
Temperature compensation	16
Temperature/pressure ratings	17

### U

Uncertainty of measurement						16
Use						. 4

#### W

Warnings				 	 		 	 	 	3
Weight										17
Workplace safety										4



71515065

## www.addresses.endress.com

